


# KLINIKK ANGEL

Aesthetics & Beauty Equipment

## 4 FREQUENCIES

Optimal energy delivery to precise area, improving treatment results and comfort experience

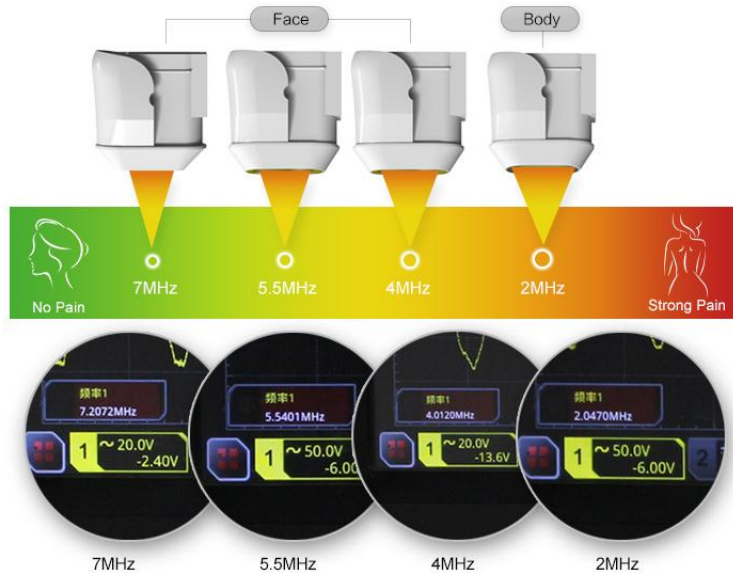

Face Lift • Skin Tightening • Body Sculpting



# LAXUS MMFU MPT

Applicable to various parts






- 1. MMFU: micro and macro together.  
MFU (micro) : less hurt, for face  
MFU (macro): deeper, for body
- MPT: ( Micro pulse technology), 416 very small dots make up the line mode for lift skin.
- Four Frequencies: different frequencies to treat different skin tissues, better result, real no pain.

## Advantages



- Different frequencies
- Latest technology
- More effective
- More comfortable





## Lasting And More Effective

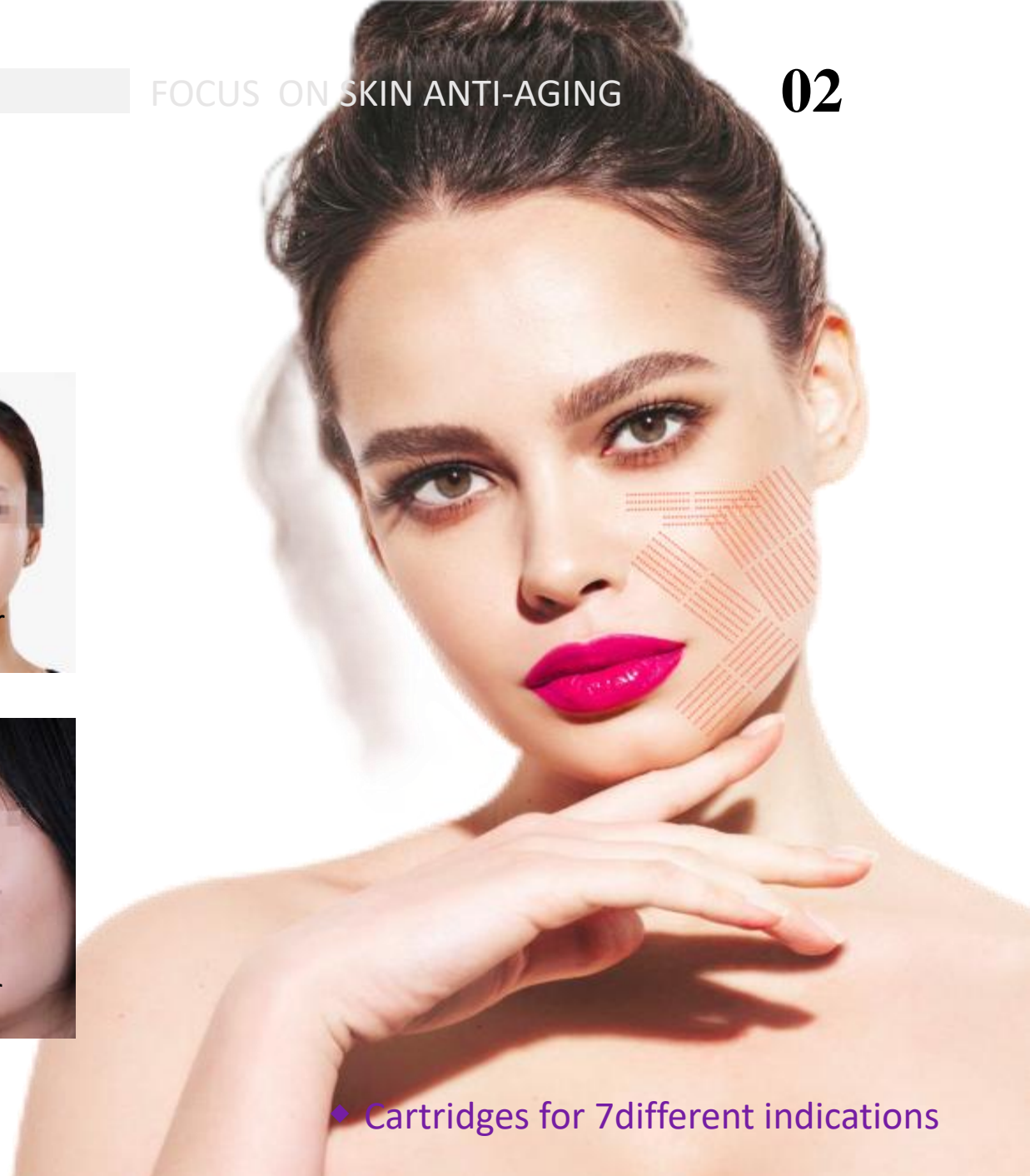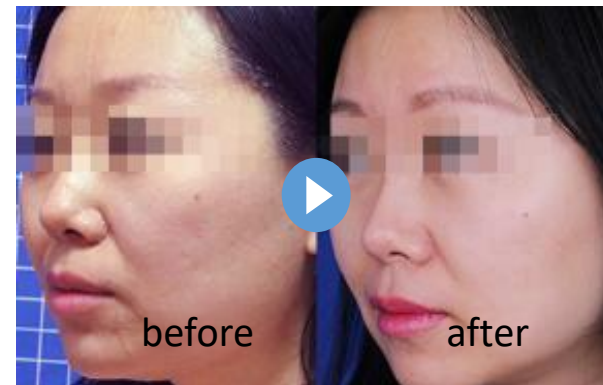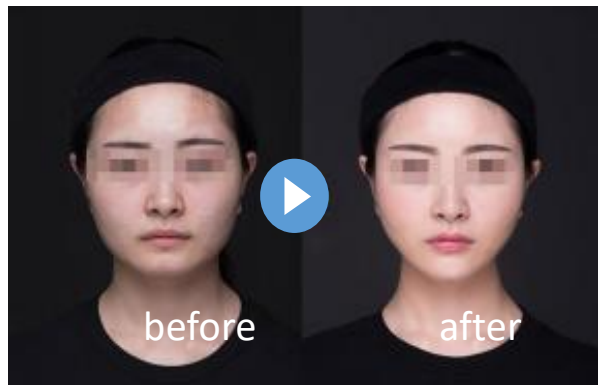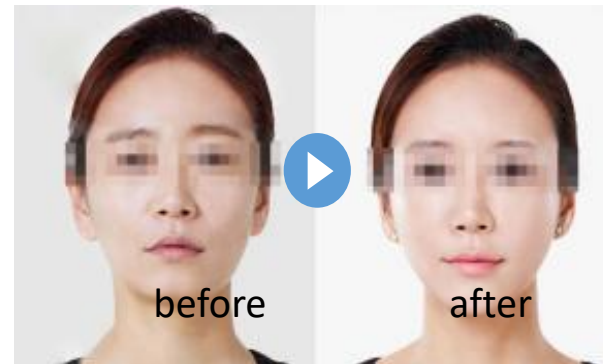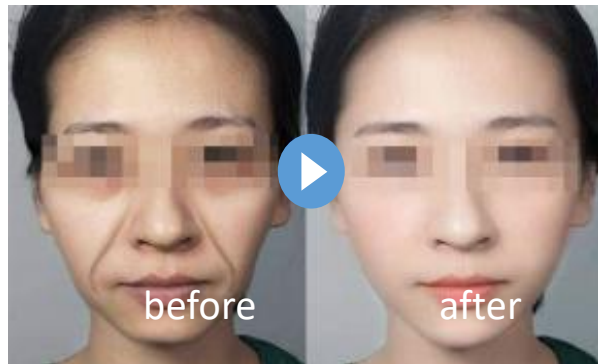

◆ Cartridges for 7 different indications

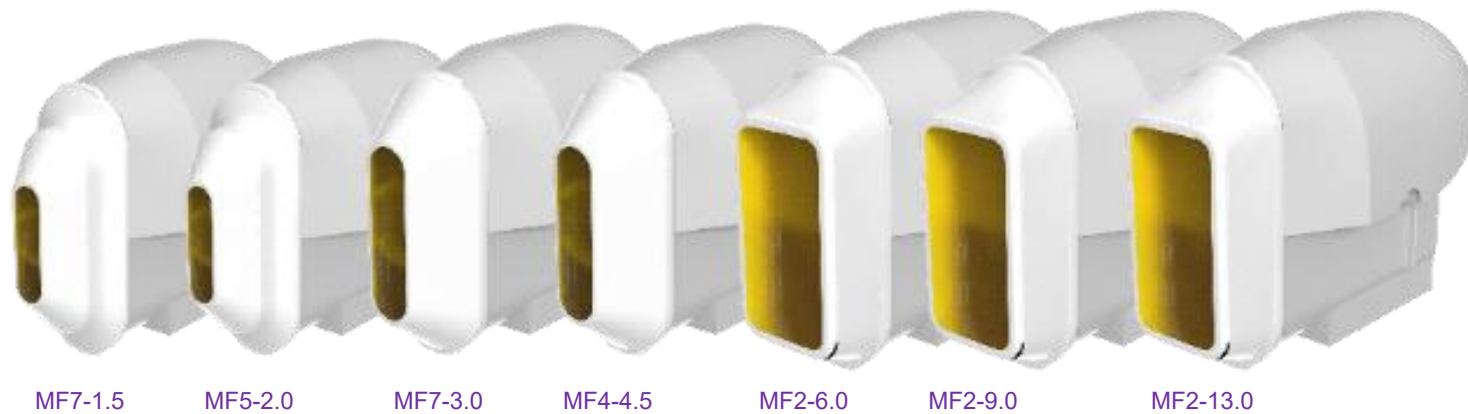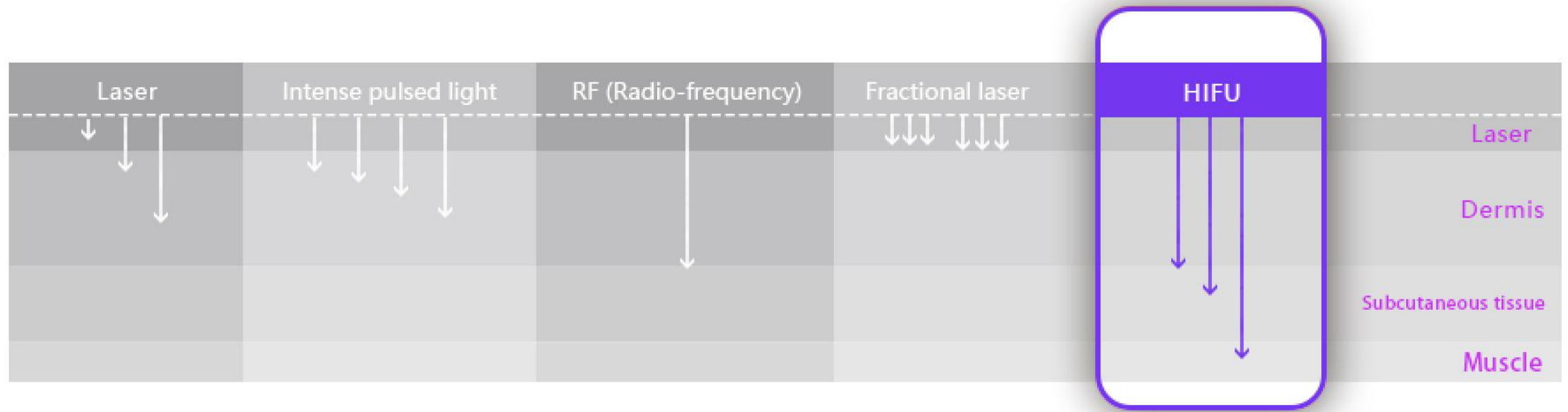



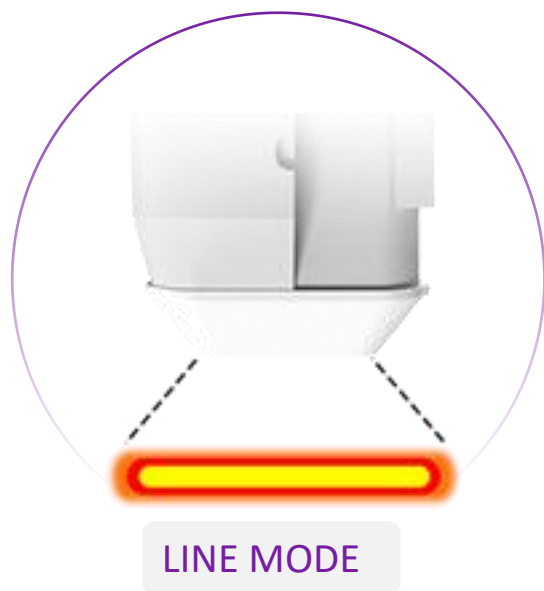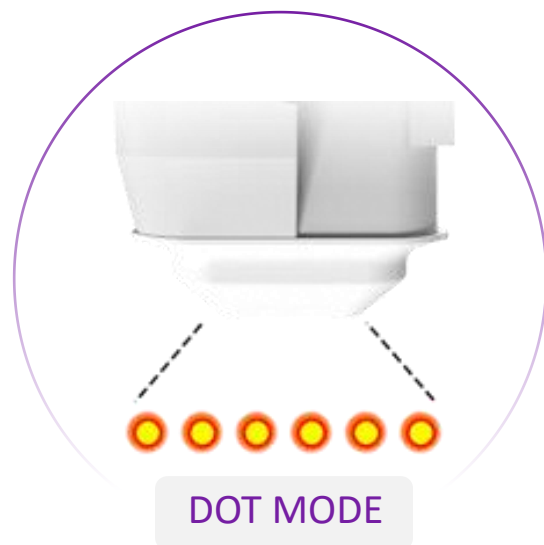

- Use Dot cartridges as a focus area with powerful energy. (Local heating) and the use of linear cartridges to evenly distribute the entire face with a wider area of energy for optimal contouring (overall heating)

Face Cartridge

MF4-4.5





MF2-13.0  
Body Cartridge




Suitable for small areas

Qs 2.0mm can also treat narrow areas like eye circles, which have a narrow tip

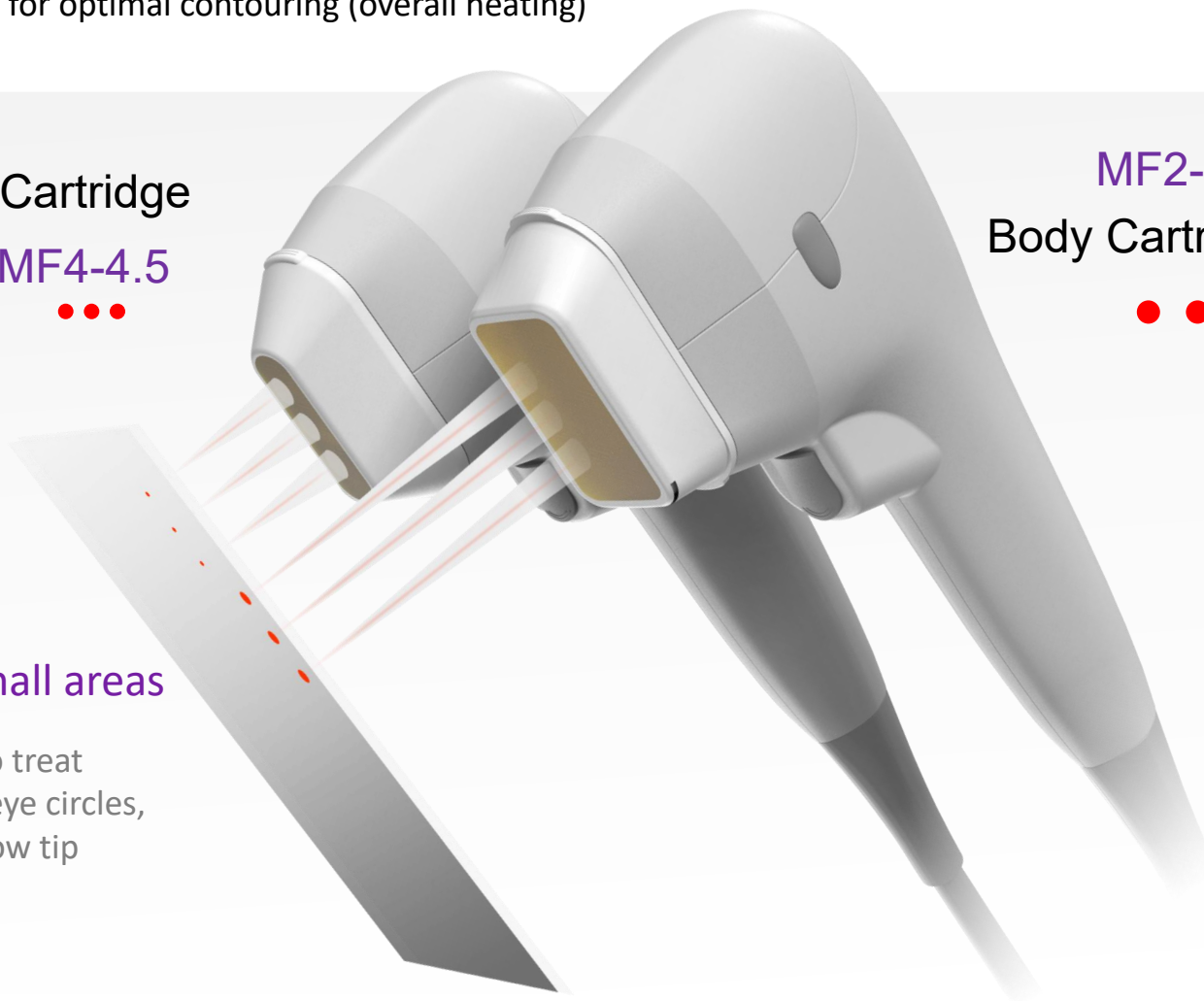


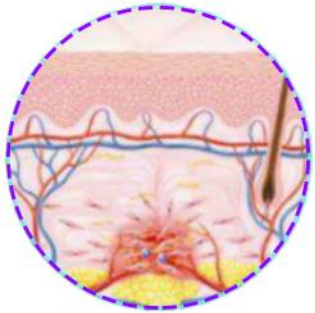

Hemostasia  
Coagulação da  
pele Migração  
celular

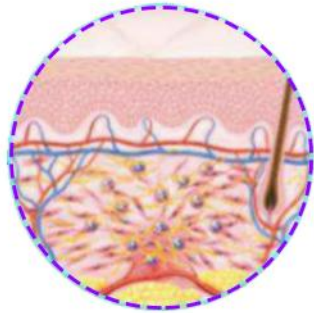

Tecido de  
granulação  
Contração de lesão  
Síntese de colágeno

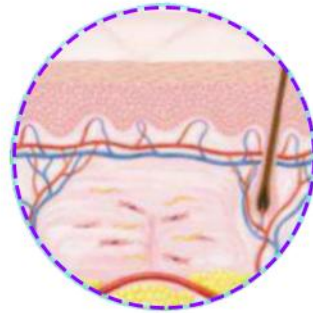

Reagrupamento do  
colágeno Aumento do tônus  
da pele Efeito lifting

### ◆ How Does the LAXUS MMFU MPT fat removal ?


#### Introduction


High Intensity Focused Ultrasound (HIFU) directly delivers heat energy to skin and subcutaneous tissue that can stimulate and renew the skin's collagen and thus consequently improving the texture and reducing sagging of the skin. It literally achieves the results of a facelift or a body lifts without any invasive surgery or injections, moreover, an added bonus of this procedure is that there is no downtime.

### ◆ How LAXUS MMFU MPT Skin Tightening works?

#### Introduction

HIFU has a smaller focus point than other. With more accurately transmits high-energy focused ultrasound energy at  $65\sim 75^{\circ}\text{C}$  to the target skin tissue layer, which results in a hot coagulation effect without harming surrounding tissues. While stimulating the proliferation of collagen and elastic fibers, it greatly improves the comfort and give you a perfect face with the skin firm, and face tightening

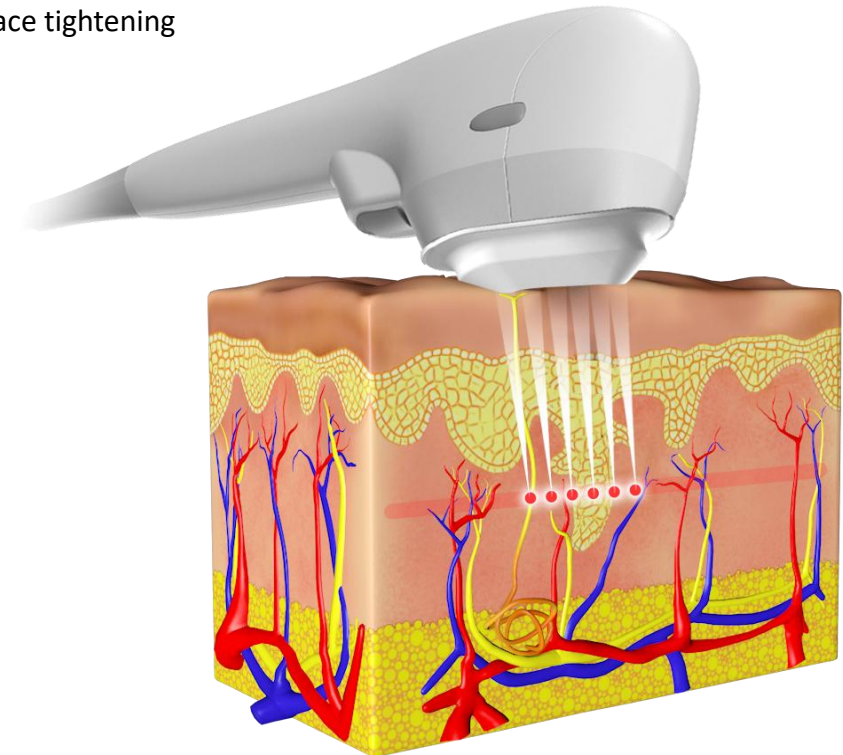

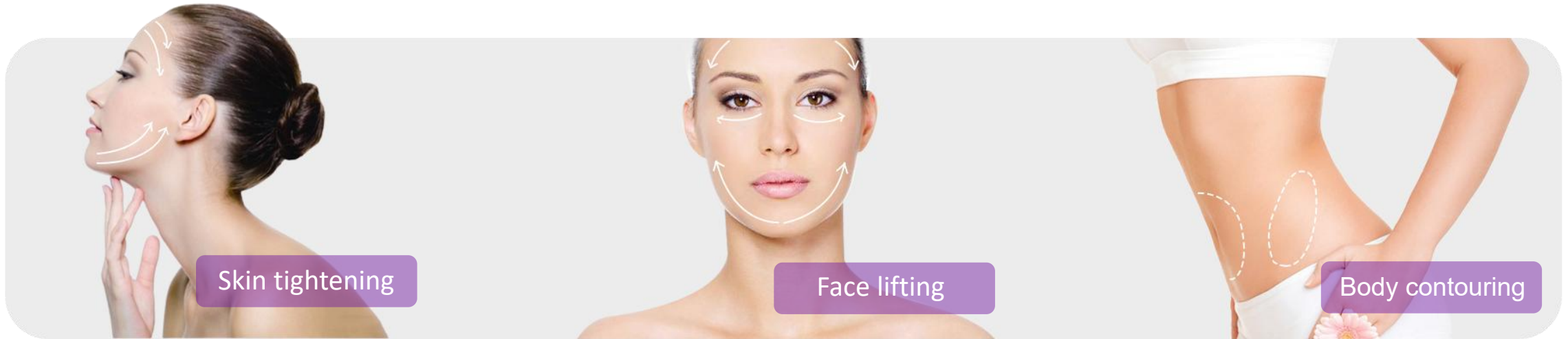


## 01 Lift & Tighten

Four multi-depth cartridges directly transfer ultrasound energy into multiple layers of the face to trigger collagen renewal for a tighter, younger complexion. The 2.0 mm cartridge can also treat both narrow and sensitive contours with a slimmer transducer and faster shot speed.



## 02 Tighten & Body Sculpting

Macro-focused cartridges at 6.0mm, 9.0mm, and 13.0mm directly transfer an equal distribution of ultrasound energy into the subcutaneous tissue of multiple body areas, inducing contractions as well as naturally disposing of stubborn fat cells for contouring purposes.

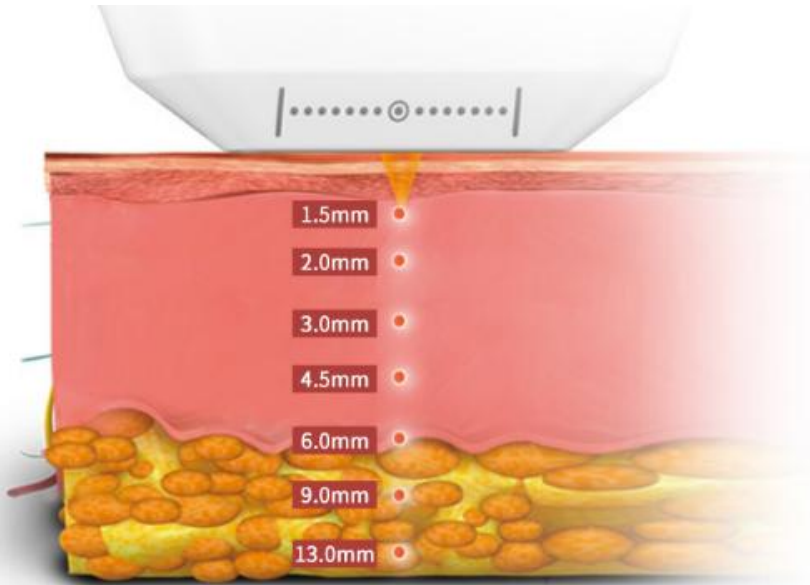


## 4 frequencies Precision treatment

Painless due to optimal energy delivery to the exact target area.

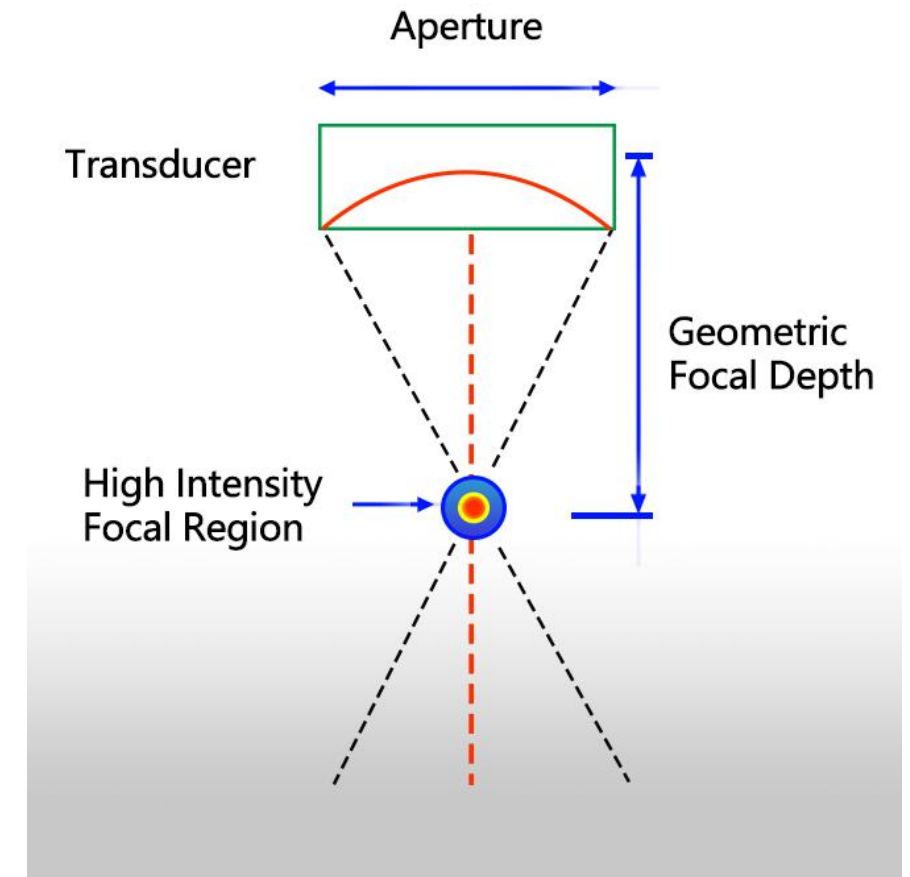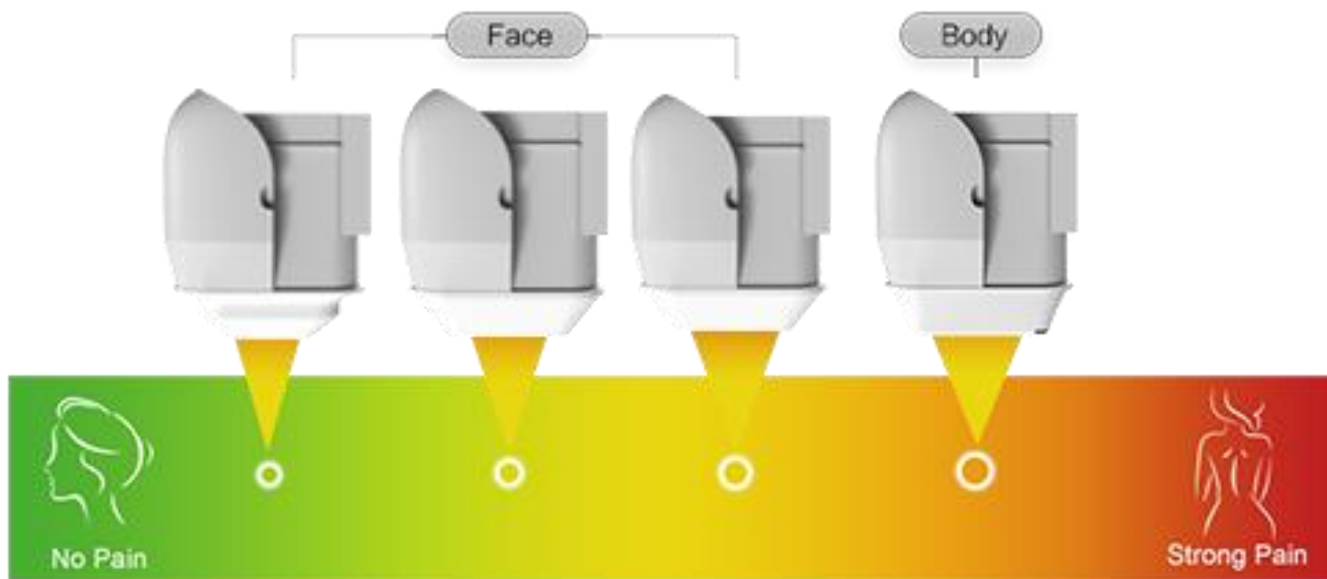



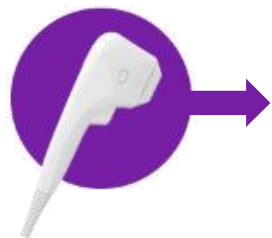


## ◆ Face




## ◆ Body



## ◆ Face

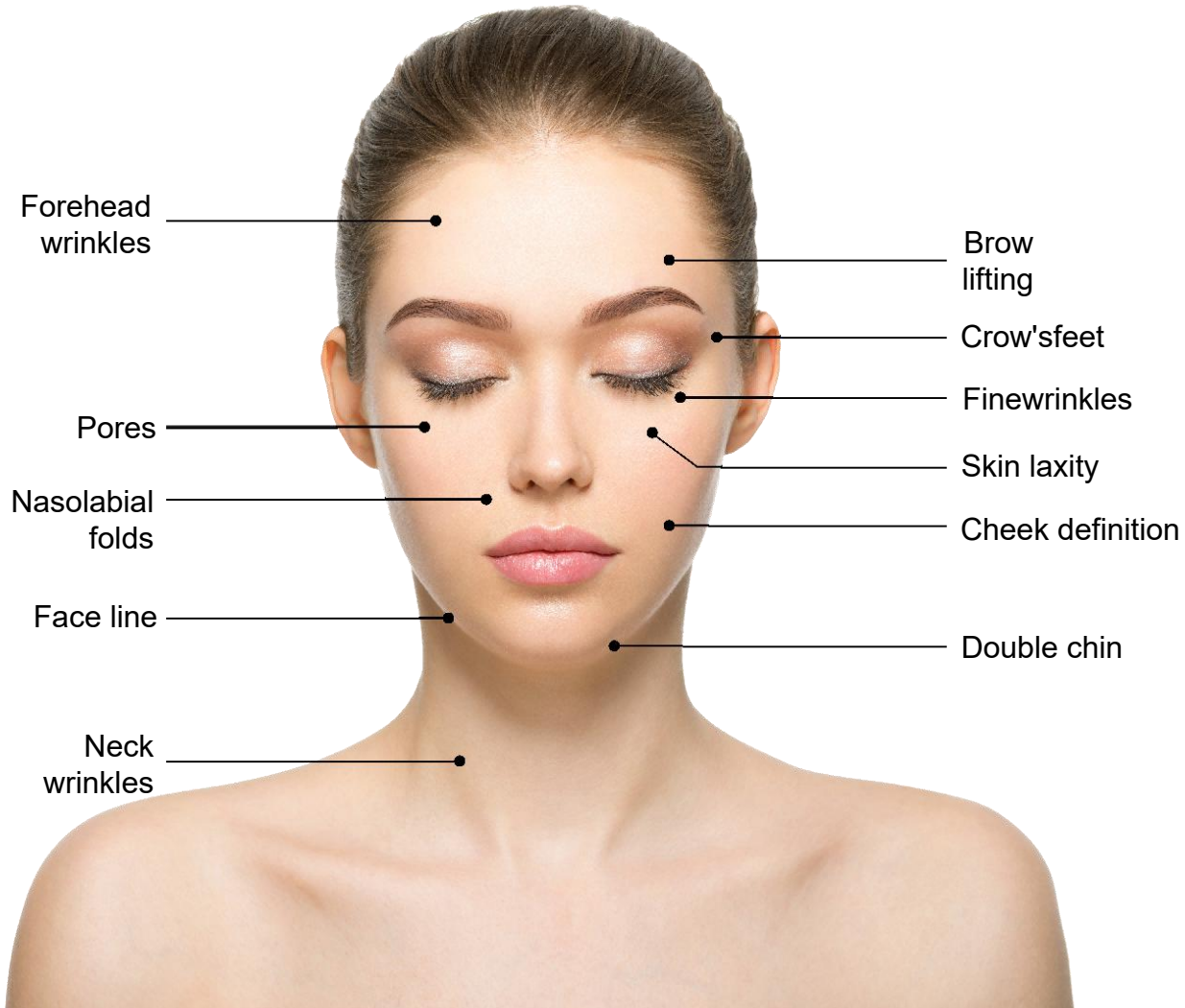

## ◆ Body

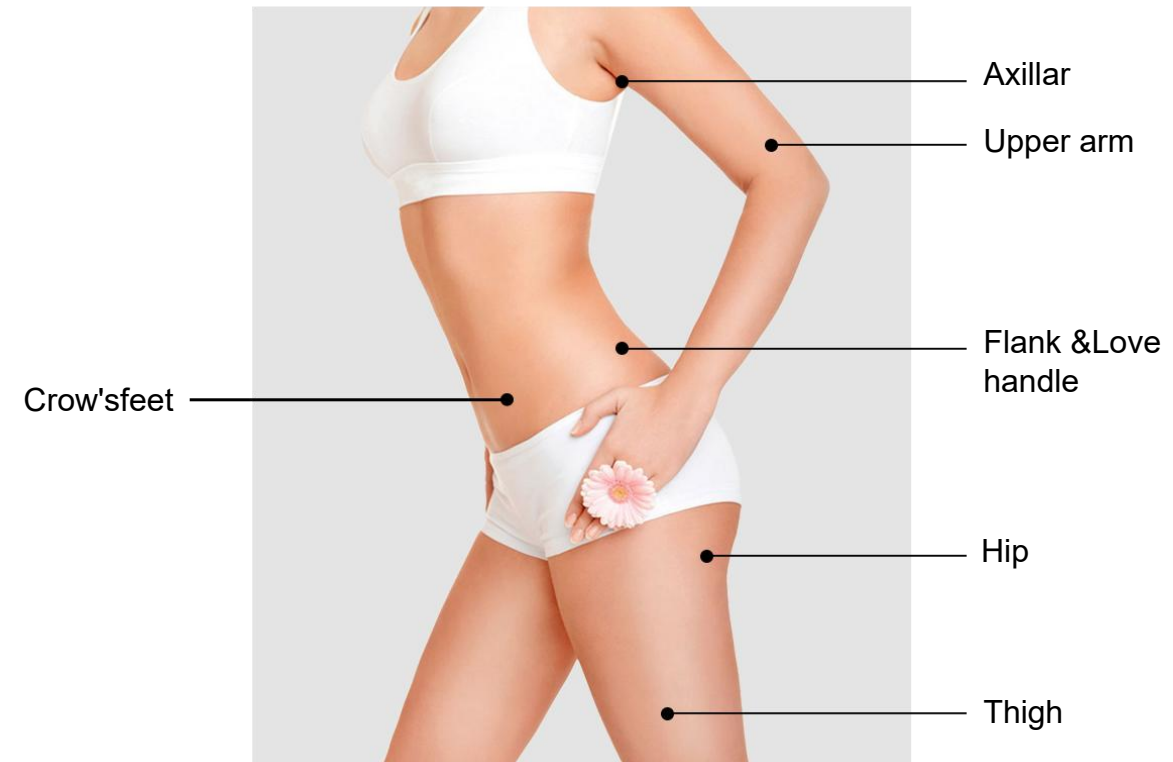

## FACE



Different treatment heads are more effective for different parts of the body



Calf



Chest



buttocks



Abdomen



Thigh



Both side of waist




Back



Arm

## BODY

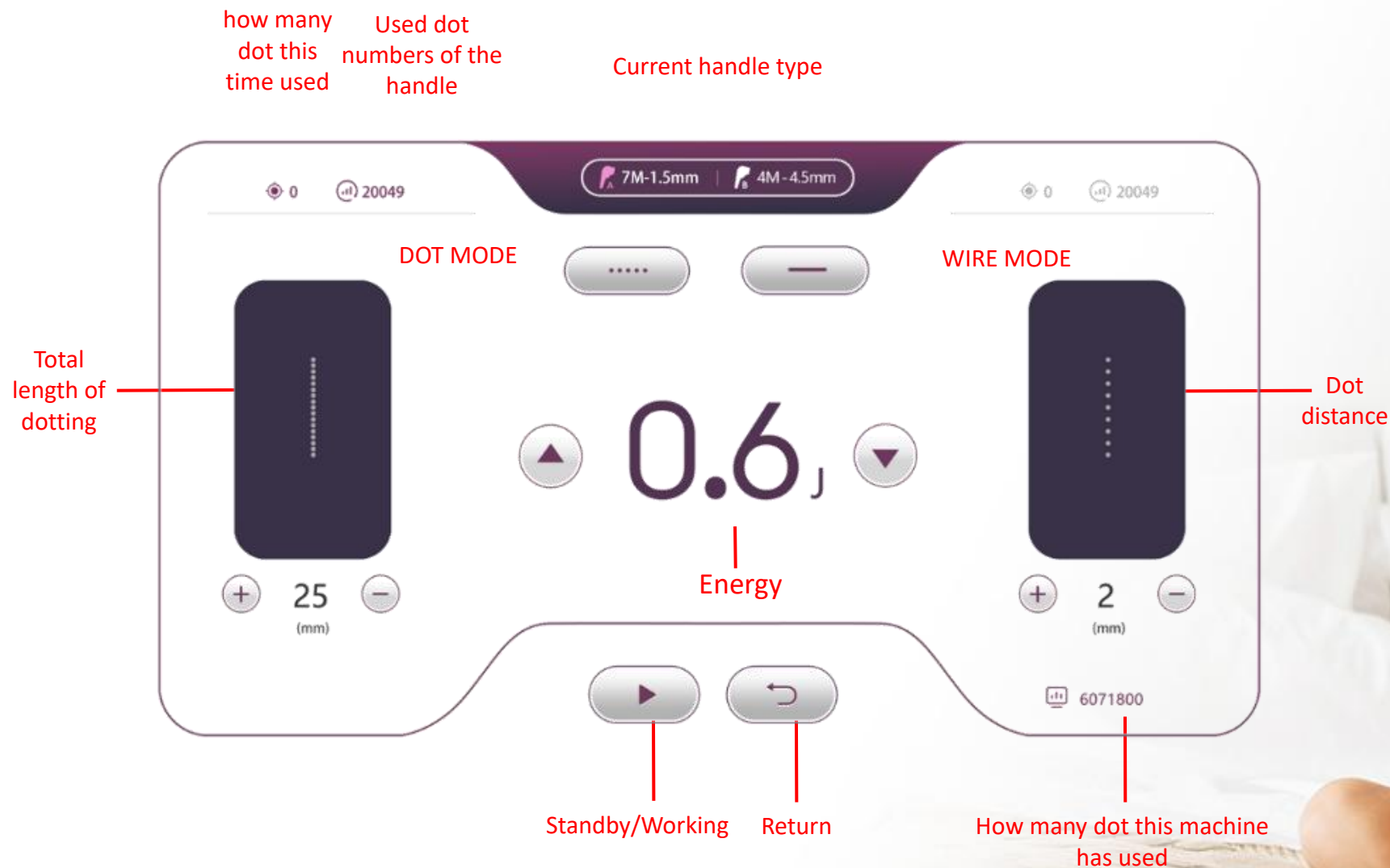



# Ancillary product display



	1.5mm [7M, 0.1-3.0J]	
	2.0mm [5M, 0.1-3.0J]	
	3.0mm [7M, 0.1-3.0J]	
	4.5mm [4M, 0.1-3.0J]	
	6.0mm [2M, 0.1-3.0J]	
	9.0mm [2M, 0.1-3.0J]	
	13.0mm [2M, 0.1-3.0J]	
shot times	30000	
spacing	1.0-2.0mm	Stepping 0.1mm
length	5-25mm	Stepping 1mm
Power supply	AC100-240V 50/60Hz	



Package box

Treatment head

Handle

# Thanks!

## CONTACT US

Klinikk Angel AS  
Aesthetics & Beauty Equipments

**Melhusvegen 505**  
**7224 Melhus**  
**Norway**

---

Tel + 47 979 78 941

Web: [www.klinikkangel.no](http://www.klinikkangel.no)

Email: [post@klinikkangel.no](mailto:post@klinikkangel.no)